


ISWIM & ISWIM 2

The iSwim is a product developed by Autolift, it is a portable pool lift that has been designed to allow users to get into a pool at any location with ease.

Features include;

- Portable – can enter and exit the pool at any location depending on the size of the coping used on the side of the pool
- Anti-Tip and Stabilizing System
- Lifting capacity 136kg
- Emergency Stop Button

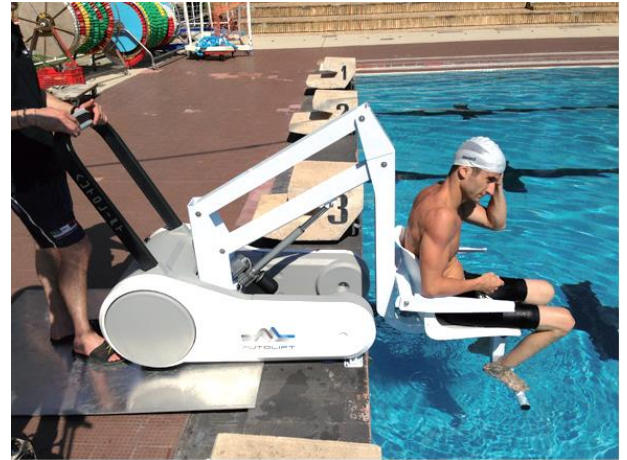

iSwim

Structure Materials – Aluminium/ Stainless Steel
 Wheels – 4 (2 with brake)
 Unit weight 140kg
 Battery – 1 x 30 Amp
 Battery Working Range 50 cycles approx.
 Battery Charger – Included
 Length – 1245mm
 Width – 640mm
 Height – 1510mm
 Rotation - Nil
 Functionality on Handrail – Battery Load indicator, up-down push buttons, emergency stop button

iSwim 2

Structure Materials – Aluminium/ Stainless Steel
 Wheels – 4 (2 with brake)
 Unit Weight 300kg
 Battery – 2 x 18 Amps
 Battery Working Range 30 cycles approx..
 Battery Charger – Included
 Length – 1245mm
 Width – 640mm
 Height – 1510mm
 Rotation – 200 degrees
 Functionality on Handrail – Battery Load indicator, up-down push buttons, emergency stop button



Evolution of the iSwim

As you can see there are two options available for the iSwim, the iSwim and the iSwim 2. The iSwim 2 was developed to allow independent swimming pool access for users. Once the lift is positioned the user can make the transfer into the pool completely independently. This is because the base of the iSwim 2 has up to 200 degrees’ rotation as well as a hand held remote, allowing the user to transfer from mobility device at one angle and enter the swimming pool at another, as suitable for the user. The rotation that the iSwim 2 allows is also beneficial for some users who may find it more comfortable to transfer into the pool side ways for safety, giving the iSwim more flexibility of use.